

## **IEA Medical Migränemaske**

### **Anwendung:**

- Cold Therapy (Kältetherapie):
  - Die Maske für maximal 2 Stunden in den Gefrierschrank/Kühlschrank legen.
  - Maximale Anwendung: 2–4 Stunden in kurzen Intervallen von 20 Minuten.
- Warm Therapy (Wärmetherapie):
  - In der Mikrowelle bei 600 Watt für etwa 30–40 Sekunden erwärmen.
  - Zwischendurch prüfen, um Verbrennungen und Überhitzung zu vermeiden.

### **Warnhinweise:**


- Überhitzung und Verbrennungen vermeiden, Maske beim Erwärmen beobachten.
- Nicht länger als empfohlene Zeit anwenden (max. 20-Minuten-Intervalle).
- Außerhalb der Reichweite von Kindern aufbewahren.
- Vor der Nutzung auf Hautreaktionen achten, bei Beschwerden sofort absetzen.

### **Kontakt-Hersteller:**

- [iris@handasen.com](mailto:iris@handasen.com)

Use as warm therapy

Alleviates: Headache, migraine

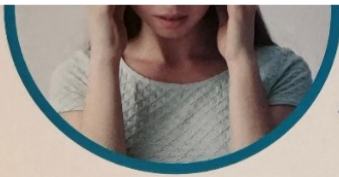

Soft material - fits perfectly head shape

Light inhibiting - Blocks light irritation

APPLICATION



**COLD THERAPY**  
Put in the freezer/refrigerator for a maximum of 2 hours; Maximum application 2-4 hours in short intervals of 20 minutes.



**WARM THERAPY**  
Microwave at 600 watts for about 30-40 seconds. Check in between to avoid burns and overheating.



COLD THERAPY

WARM THERAPY



IEA International Trading GmbH  
Grünwalder Straße 61  
81547 München

info@iea-medical.com

+49 (0) 89-45 20 503-0

REF 20240126


LOT 1B260124




JIANGXI HANDASEN INDUSTRY CO., LTD.  
Shaxi Avenue, Zhennan New District, Shaxi Town,  
Xinzhou District, Shangrao City, Jiangxi, China

EC REP

Riomavix S.L., Calle de Almansa 55, 10,  
Madrid 28039 Spain

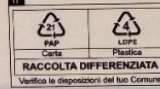

25.01.2026

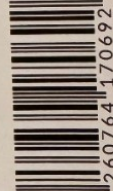

4 260764 170692




1 MASKE



# MIGRÄNE MASKE

Migraine mask  
Masque anti-migraineux  
Maschera per l'emicrania  
Máscara de migraña

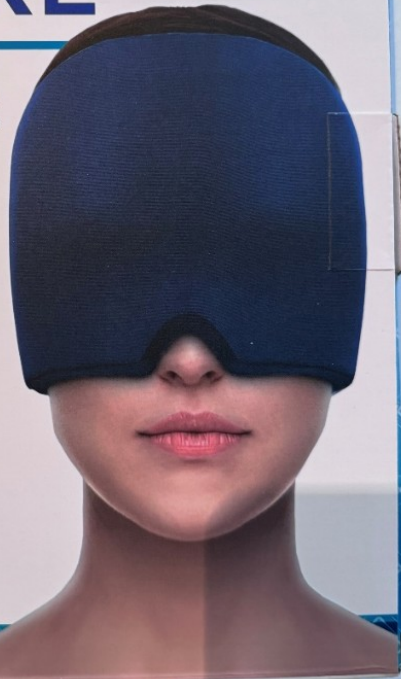

EASY TO USE



REUSABLE



HIGH COMFORT



UNIVERSAL FIT  
/ADAPTS TO THE SHAPE OF  
THE HEAD INDIVIDUALLY

## APPLICATION AREAS



Use as cold therapy



Use as warm therapy



Alleviates: Headache, migraine



Promotes: Circulation (Warm)  
Firms the skin (Cold)



Soft material - fits perfectly head shape



Light inhibiting - Blocks light irritation

## APPLICATION



### COLD THERAPY

Put in the freezer/refrigerator for a maximum of 2 hours. Maximum application 2-4 hours in short intervals of 20 minutes.



### WARM THERAPY

Microwave at 600 watts for about 30-40 seconds. Check in between to avoid burns and overheating.



COLD THERAPY

WARM THERAPY



IEA International Trading GmbH  
Grünwalder Straße 61  
81547 München

info@lea-medical.com

+49 (0) 89-45 20 503-0

REF 20240126


LOT 1B260124



JIANGXI HANDASEN INDUSTRY CO., LTD.  
Shaxi Avenue, Zhennan New District, Shaxi Town,  
Xinzhou District, Shangrao City, Jiangxi, China



Riomarx S.L. Calle de Almansa 55, 10,  
Madrid 28039 Spain

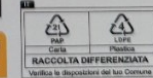

4 260764 170692

25.01.2026